

Operating Instructions ■

Perth® II Recliner

September 2011

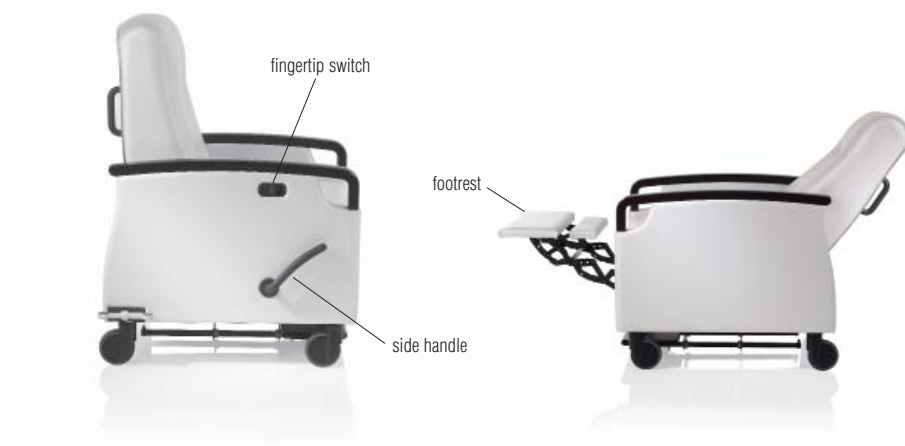


Figure 1

Recline Operation

To engage the back recline, pull back on the fingertip switch while leaning back until you reach the desired degree of recline. When the desired recline position is met, release your fingers to stop the recline motion. If the recline needs to be adjusted, pull back on the fingertip switch and lean backward or forward to the desired position, then release. To move the chair back into the sitting position, pull back on the fingertip switch, lean forward and the back will follow you to an upright position.

Footrest Operation

Pull back on the side handle while in the sitting position to raise the footrest. Push forward on side handle while in the sitting position to close the footrest (Figure 1). Do not step on footrest to close.

Note: Stepping on the footrest to close will shear off the activator tube that is connected to the side handle. It will torque the metal inside because of excessive force and cause it to break off with the handle. The user, in the sitting position of the recliner, should be using only the handle to open/close the chair.