

Kindera Recliner Order Check List:

1 Model Number

2 Upholstery Selection

Please supply the following:

Fabric vendor/Pattern name/Pattern color

[View all JSI Standard and Graded-In Textiles](#)

- Indicate direction of fabric - Railroaded or Up the Roll
- Please provide illustration, if available

All COM orders must include the following:

- Vendor's name, pattern name and color
- Sample of the fabric attached to the order
- Model numbers of items for fabric application
- Special application instructions
- Customer and purchase order number


Ship prepaid COM material to:

JSI
Attn: COM Department
201 East Martin Street
Receiving B, Dock 71
Orleans, IN 47452

Application of fabric:

Railroaded

Up the Roll

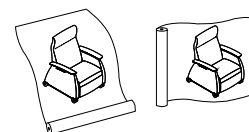


All COM orders will ship at standard lead time from receipt of COM except for orders requiring special custom matched finishes.

Upholstery Application

Kindera offers the ability to be specified with multiple upholsteries to create a unique look. Both upholstery styles are presented in the price list and eTools. Models not ending in a _M will be upholstered in the single upholstery for the complete unit.

Models ending with a _M designate multiple upholstery zones will be specified (ex: 43SPFA_M) To calculate pricing, add base price + zone 1 upcharge + zone 2 upcharge + zone 3 upcharge.

When specifying multiple upholstery, make sure to identify specific yardage per zone.

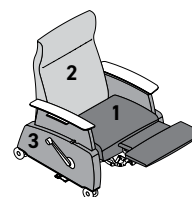

Zone 1 = Seat & Footrest

Zone 2 = Back

Zone 3 = Arms with Back Kickplate

3 Moisture Barrier

Offers further protection from fluids:

NMB No Moisture Barrier

MB1 Moisture Barrier Seat Only - upcharge applies

MB2 Moisture Barrier Seat & Back - upcharge applies

4 Transfer Arm

Kindera features the patent pending Trans4rm™ system in both the dual transfer and single transfer arm models.

Dual transfer recliners allow either arm to retract into the recliner frame, creating a seamless way to transfer patients without increasing the footprint. Single transfer models provide one fixed arm and one transfer arm, enabling specification of left or right for transfer. Models using the Trans4rm™ system smoothly retract their arm into the frame, creating a smooth surface for safe patient transfers.

5 Arm Caps

Arm Caps are available in three different materials:

Polyurethane - standard

SGU Stone Grey

3D Laminate - upcharge applies

BEL Belair
FLO Florence Walnut
PIL Pilsner
PIN Pinnacle Walnut

Solid Surface - upcharge applies

DWH Designer White



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Recliner Order Check List (continued):

6 Footrest Handle Location

Footrest handle is standard on all models with placement on the right as seated by default.

RH Footrest handle on right
LH Footrest handle on left
LRH Footrest handle on both left and right – upcharge applies

7 Push Handles

To aide in repositioning of chair optional push handles can be specified as a set of two on the recliner or ordered individually as a replacement part.
Note: Recliner is not intended for patient transport.

NPH No push handles
YPH Yes push handles – upcharge applies

8 Brake/Casters

Kindera is equipped with (4) 4" heavy-duty, non-marking, dual-wheel swivel casters. Standard Caster/Brake configuration includes (2) front individual locking and (2) back non-locking casters.

All models offer central lock as an upgrade. This option features dual-sided levers positioned at arm height, to easily lock or release the casters. In this configuration, the left front and right rear casters (as seated) are locking, while remaining two casters are non-locking.

2CB 2 Locking Caster Brake
CNB Central Caster Lock Lever – upcharge applies

9 Heat & Massage

Optional Heated Seat & Back and Massage enhance patient comfort.

On fixed arm units, the control panel is located on the right arm; on single transfer models, it is located on the non-transfer arm. This option is not available on dual transfer models. The power cord exits the bottom of the chair, and when Heat, Massage, and Power are selected together, a single combined cord is used. Heat is regulated to not exceed 102°F. Fixed arm units include three heat settings, with optional massage for enhanced long-term comfort.

NHSB No Heat or Massage
YHSB Yes Heated (Seat & Back) – upcharge applies
YHMSB Yes Heat & Massage (Seat & Back) – upcharge applies

10 Foley Hook

Foley hook available with single footrest handle; mounts opposite handle. Not available with folding table. Hook is always black and features magnetic closure.

NFH No Foley Hook
YFH Foley Hook – upcharge applies



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Recliner Order Check List (continued):

11 Power in Arm

A USB A+C black power port with integrated spill shield is available for patient convenience.

In fixed arm units, the control panel is always on the right arm. In single transfer models the control panel is always on the non transfer arm. This option is not offered on dual transfer arm models. The cord extends from the bottom of the chair, if heated seat and power are selected on same chair a combined power cord is used.

Fixed Arm Models

NPU No Power
 PUL Power in Left Arm (as seated) – upcharge applies
 PUR Power in Right Arm (as seated) – upcharge applies

Single Transfer Arm Models

NPU No Power
 YPU Power in Fixed Arm – upcharge applies

12 Foot Tray

A roll-out black foot tray with durable KYDEX® thermoplastic surface delivers added support when needed, then tucks away to keep the space clear. Integrated foot pull allows hands-free extension.

width = 14 3/8"
 maximum extension = 13 3/4"

NFY No Foot Tray
 YFY Yes Foot Tray – upcharge applies

13 Folding Table

Designed for convenient one-handed use, the folding table locks securely in place, keeping your essentials easily accessible.

Note: Not available on the side with the footrest handle.

3D Laminate – upcharge applies

BEL Belair
 FLO Florence Walnut
 PIL Pilsner
 PIN Pinnacle Walnut

Solid Surface – upcharge applies

DWH Designer White
 HST Hot Stone
 ONX Black Onyx Mirage
 SOG Soothing Grey



CORPORATE OFFICE
 225 CLAY STREET
 JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
 812.482.3204 **OFFICE**
 812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Recliner Specifications

1. Trans4rm™ System

Allows arm to retract into the recliner's frame, creating a seamless way to transfer patients without increasing the recliner's footprint. Available as single or dual.



2. Dual Back Recline Levers

Accessible, elongated levers below the arm caps easily activate the back's infinite recline.



3. Trendelenburg (Dual-Sided)

Foot pedal at the base quickly repositions the reclined chair to elevate the patient's feet above their head.



4. Heat & Massage

Fixed arms offer three heat settings and optional massage for added comfort over time.



5. Foot Tray

A roll-out black foot tray with durable KYDEX® thermoplastic surface delivers added support when needed, then tucks away to keep the space clear. Integrated foot pull allows hands-free extension.



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Recliner Specifications

Reclined Positions



Bariatric

infinite recline 60°

footprint for full back/footrest extended is 64.8"



Patient

infinite recline - 25°

footprint for full back/footrest extended to lay flat is 74.5"



Trendelenburg

available in standard size only

activated by the foot pedal along the base

Clean Out

Ensures simple and thorough cleaning around the seat.



Recline Motion/Weight Capacity/ Cleanable Specs

- Each Kindera model has infinite recline from sitting position to 25° for patient model and 60° for bariatric models.
- Kindera is thoughtfully designed to accommodate the regimens of clinical environments with durable arm caps, bumper guard, recline lever, footrest handle and push handles that wipe down quickly and clean outs to ensure ease of cleaning.
- Built for demanding healthcare facilities, Kindera exceeds rigorous testing.
ANSI/BIFMA X5.41 (400 lb)-2021 - Large Occupant Public and Lounge Seating Tests¹
ANSI/BIFMA X5.41 (600 lb)-2021 - Large Occupant Public and Lounge Seating Tests²

Weight Capacities are as follows:

500 lb. Patient Recliner
500 lb. Glider
750 lb. Bariatric Recliner
80 lb. Breakaway ottoman



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Glider Order Check List:

1 Model Number

2 Upholstery Selection

Please supply the following:

Fabric vendor/Pattern name/Pattern color

[View all JSI Standard and Graded-In Textiles](#)

- Indicate direction of fabric - Railroaded or Up the Roll
- Please provide illustration, if available

All COM orders must include the following:

- Vendor's name, pattern name and color
- Sample of the fabric attached to the order
- Model numbers of items for fabric application
- Special application instructions
- Customer and purchase order number


Ship prepaid COM material to:

JSI
Attn: COM Department
201 East Martin Street
Receiving B, Dock 71
Orleans, IN 47452

Application of fabric:

Railroaded

Up the Roll

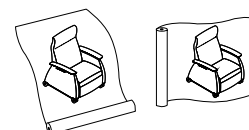


All COM orders will ship at standard lead time from receipt of COM except for orders requiring special custom matched finishes.

Upholstery Application

Kindera offers the ability to be specified with multiple upholsteries to create a unique look. Both upholstery styles are presented in the price list and eTools. Models not ending in a _M will be upholstered in the single upholstery for the complete unit.

Models ending with a _M designate multiple upholstery zones will be specified (ex: 43SPGFA_M) To calculate pricing, add base price + zone 1 upcharge + zone 2 upcharge + zone 3 upcharge.

When specifying multiple upholstery, make sure to identify specific yardage per zone.

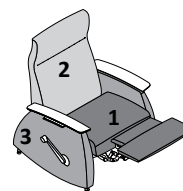

Zone 1 = Seat & Footrest

Zone 2 = Back

Zone 3 = Arms with Back Kickplate

3 Moisture Barrier

Offers further protection from fluids:

NMB No Moisture Barrier

MB1 Moisture Barrier Seat Only - upcharge applies

MB2 Moisture Barrier Seat & Back - upcharge applies

4 Arm Caps

Arm Caps are available in three different materials:

Polyurethane - standard

SGU Stone Grey

3D Laminate - upcharge applies

BEL Belair
FLO Florence Walnut
PIL Pilsner
PIN Pinnacle Walnut

Solid Surface - upcharge applies

DWH Designer White

5 Glide Lock Handle

Allows users to transition between gentle motion and a stable, stationary position. Located on the left as seated.



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Glider Features

1. Dual Back Recline Levers

Accessible, elongated levers below the arm caps easily activate the back's infinite recline.



2. Glide Lock Handle

Provides a stationary position when engaged.



3. Adjustable Leveling Glides



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI

Kindera Replacement Parts

Easy Replacement Parts

Designed for healthcare's challenges, Kindera recliner ensures easy replacement parts if damaged

Replaceable high-use parts extend the product's lifespan

This includes: seat, mid-ottoman, footrest, arm caps, individual locking casters, push handles, folding table top, and leveling glides, which are detailed in the price list.



CORPORATE OFFICE
225 CLAY STREET
JASPER, INDIANA 47546



800.457.4511 **TOLL FREE**
812.482.3204 **OFFICE**
812.482.1548 **FAX**



JSIFURNITURE.COM



CONTACT JSI